


**Kwik-Skan®**

The Kwik-Skan monitors feature a unique oval shape that optimizes contact with vessels in the forehead for accurate monitoring. Temperature is indicated by bar and graphical scales in both Fahrenheit and Celsius.

Temperature strips record temperatures by averaging the temperature of the tissue under the strip. Therefore the size and shape of the temperature strip can significantly affect accuracy. Kwik-Skan's vertically oriented oval design mimics the main vessels of the forehead for better accuracy.

- Unique Oval shape
- Monitors temperatures from 92° to 106°F (34° to 41°C)
- Bar and Numeric scales simplify reading temperatures
- Available in boxes of individual envelopes or in a bulk dispenser box
- Adhesive strips on boxes for easy mounting to walls, cabinets or carts
- Can be used on the forehead or anywhere temperature monitoring is desired

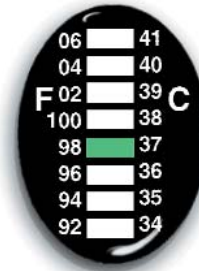

(Actual size)



Order #	Description	Unit of Sale	Box Count
5030	Kwik-Skan Oval with Dual Scale - Individual Packages 92° - 106°F / 34° - 41°C	box	100
5130	Kwik-Skan Oval with Dual Scale - Dispenser Box 92° - 106°F / 34° - 41°C	box	100